

Using A Flexhaler®

Many medications are available as inhaled treatments. Inhaled methods deliver medication directly to the airway, which is helpful for lung diseases. The patient and health care provider can choose from a variety of delivery systems for inhaling medication.

A Flexhaler® consists of a plastic device containing dry powdered medication. A Flexhaler® is breath activated. This means when you inhale, the Flexhaler® automatically releases the medication. When inhaled correctly, the medication has a better chance to reach the small airways. This increases the medication's effectiveness. The portable size, efficiency and convenience make the Flexhaler® a desirable method for inhalation treatment.

The following steps outline the best method to use a Flexhaler®.


Before Using a New Flexhaler®

1. Remove the cap from the Flexhaler® by unscrewing the cap.
2. Hold the Flexhaler® with the mouthpiece up.
3. Turn the brown bottom of the Flexhaler® all the way to the right and back to the left. You will hear it click.
4. **Repeat steps 2-3.** Only follow steps 1-4 the first time you use your new Flexhaler®.



Taking an Inhaled Treatment

1. Remove the cap from the Flexhaler® by unscrewing the cap.
2. Hold the Flexhaler® with the mouthpiece up. This will load the medication correctly.
3. Turn the bottom all the way to the right and back to the left. You will hear it click. This will load the medication.
4. Hold the Flexhaler® away from your mouth and gently breathe out.
5. Seal your lips around the mouthpiece.
6. **Inhale rapidly and deeply.** Continue to take a full, deep breath.
7. **Resume normal breathing.**
8. **Repeat steps 1-7 when more than one puff is prescribed.**

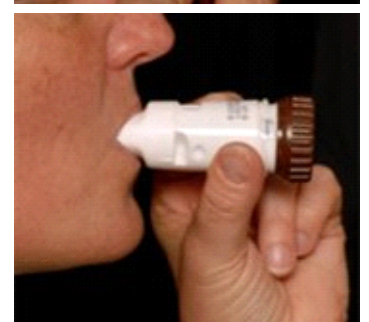

9. Rinse your mouth and spit out the water after inhaling the medication.

10. Keep the Flexhaler® cap on when not in use. This will keep the Flexhaler® clean and dry.

Identifying When the Flexhaler® is Empty:

The Flexhaler® has a dose counter window below the mouthpiece. It will let you know approximately how many doses are left. As you use the Flexhaler® the numbers will count down by 10. The Flexhaler® will count down to 0 in the center of the window when the Flexhaler® is empty. Plan to get a new Flexhaler® before you see the 0 in the window.



Cleaning the Flexhaler®:

1. Wipe the mouthpiece with a clean cloth if it appears dirty or at least once a week. Do not put the Flexhaler® in water.

Visit our website for more information about support groups, clinical trials and lifestyle information.

PTE.164 © Copyright 2007, 2009, 2014

NOTE: This information is provided to you as an educational service of National Jewish Health. It is not meant to be a substitute for consulting with your own physician.

National Jewish Health is the leading respiratory hospital in the nation. Founded 123 years ago as a nonprofit hospital, National Jewish Health today is the only facility in the world dedicated exclusively to groundbreaking medical research and treatment of patients with respiratory, cardiac, immune and related disorders. Patients and families come to National Jewish Health from around the world to receive cutting-edge, comprehensive, coordinated care. To learn more, visit njhealth.org.